

clinell®



UNIVERSAL SPRAY

Spray for surface disinfection and
cleaning of non invasive medical devices



UNIVERSAL SPRAY



CLINICAL ADVANTAGES

- **Multi-purpose**
Acts as a disinfectant and detergent, replacing the need for different types of surface cleaners.
- **The most effective universal formula spray on the market**
A mix of biocides with different mechanisms of action prevents bacterial resistance and superbug formation.
- **No damage to surfaces**
Ideal for use on mattresses and safe to use on rubber, plastics and metal.
- **Effective from 10 seconds**
- **Cost effective**
Save money by replacing several different types of products.
- **Non-flammable, no gas propellants**
- **Proven to be bactericidal, virucidal, fungicidal and tuberculocidal**
Effective against MRSA, Acinetobacter, VRE, Klebsiella, Pseudomonas, TB, Norovirus, Hepatitis B & C and many more.*

KILLS AT LEAST
99.999%*
OF GERMS AND
IS EFFECTIVE FROM
10 SECONDS.

* According to European standard tests EN1275, EN1276, EN14340, EN14561, EN14562, EN14563.

ANTIMICROBIAL ACTIVITY

Powerful, patented biocidal formulation, based on a synergistic mix of quaternary ammonium compounds and a polymeric biguanide each with a completely different mechanism of action. This ensures the widest spectrum of activity whilst mediating against resistance.

Gram-positive bacteria

Including: Staphylococcus aureus and MRSA, Streptococci, Enterococci and VRE.

Gram-negative bacteria

Including: Pseudomonas, E.Coli, Acinetobacter, Klebsiella and Salmonella.

Non-enveloped viruses

Including: Norovirus and Rotavirus.

Mycobacteria

Including: Mycobacterium terrae and Mycobacterium smegmatis.

Enveloped viruses

Including: Herpes simplex 1&2, Influenza including H1N1 & H5N1, Hepatitis B & C and HIV.

Fungi

Including: Aspergillus, Candida and Penicillium.

RANGE OF USES


Suitable for use on a wide range of medical equipment and medical surfaces such as:

Hospital beds and mattresses, wheelchairs, hoists, spinal boards, head blocks, ECG equipment and incubators.


USE DISINFECTANTS SAFELY. ALWAYS READ THE LABEL AND PRODUCT INFORMATION BEFORE USE.

ALWAYS FOLLOW MEDICAL EQUIPMENT MANUFACTURERS CLEANING PROCEDURES AND GUIDELINES.


HOW TO CLEAN SURFACES EFFECTIVELY

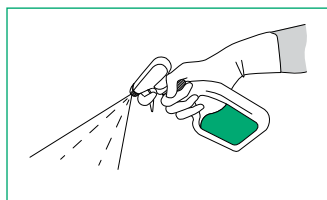

Wearing the recommended PPE, spray onto the area you wish to disinfect.

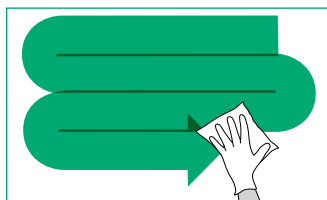

Working from clean to dirty, wipe in an 'S' shaped pattern, taking care not to go over the same area twice.

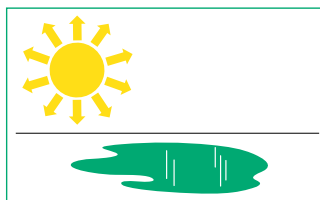

Let the surface air dry.

CE 0473

CE certified medical devices

The 500ml Universal Spray is a Class IIA medical device.

PRODUCT DESCRIPTION

Clinell Disinfectant Spray

VOLUME

500ml

NHS CODE

MFB129

ORDER CODE

CDS500